

TRANSFORMING CELL ISOLATION

EasySep™


The Easy Choice for Simple Cell Isolation



Presenter: Dr. Mausumi Ghosh

The Challenge of Cell Isolation


- Lengthy, **multi-step protocols** slow down experiments
- **High consumable** use increases lab costs
- Mechanical stress **lowers cell viability**
- User-dependent variability affects **reproducibility**

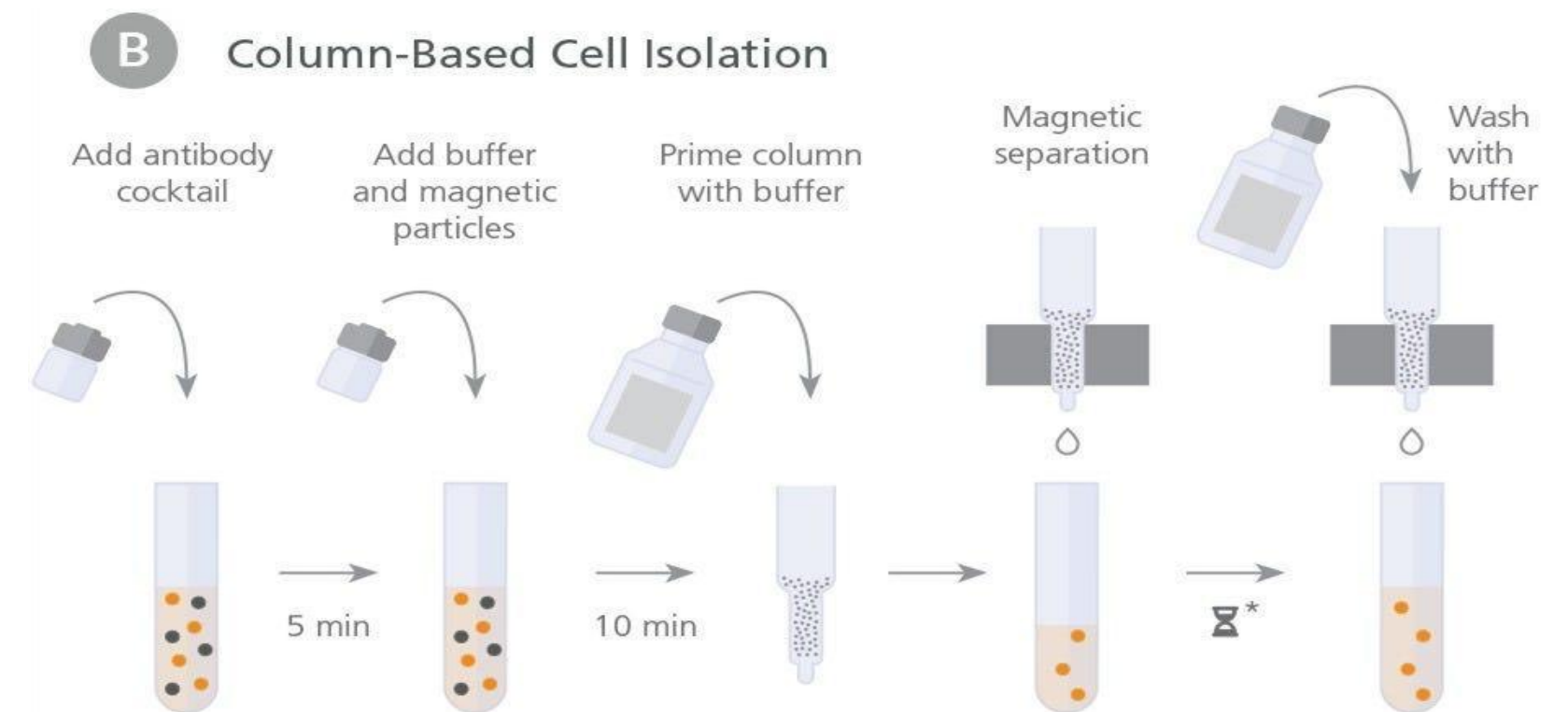

The EasySep™ Workflow

An Immunomagnetic System for Efficient Isolation

- EasySep™ is **immunomagnetic cell isolation system**.
- It uses **magnetic nanoparticles** and **antibody complexes** to separate target cells directly in the sample tube, **no columns or centrifugation**.
- Delivers **high purity, high yield**, and takes **less than 20 minutes**.

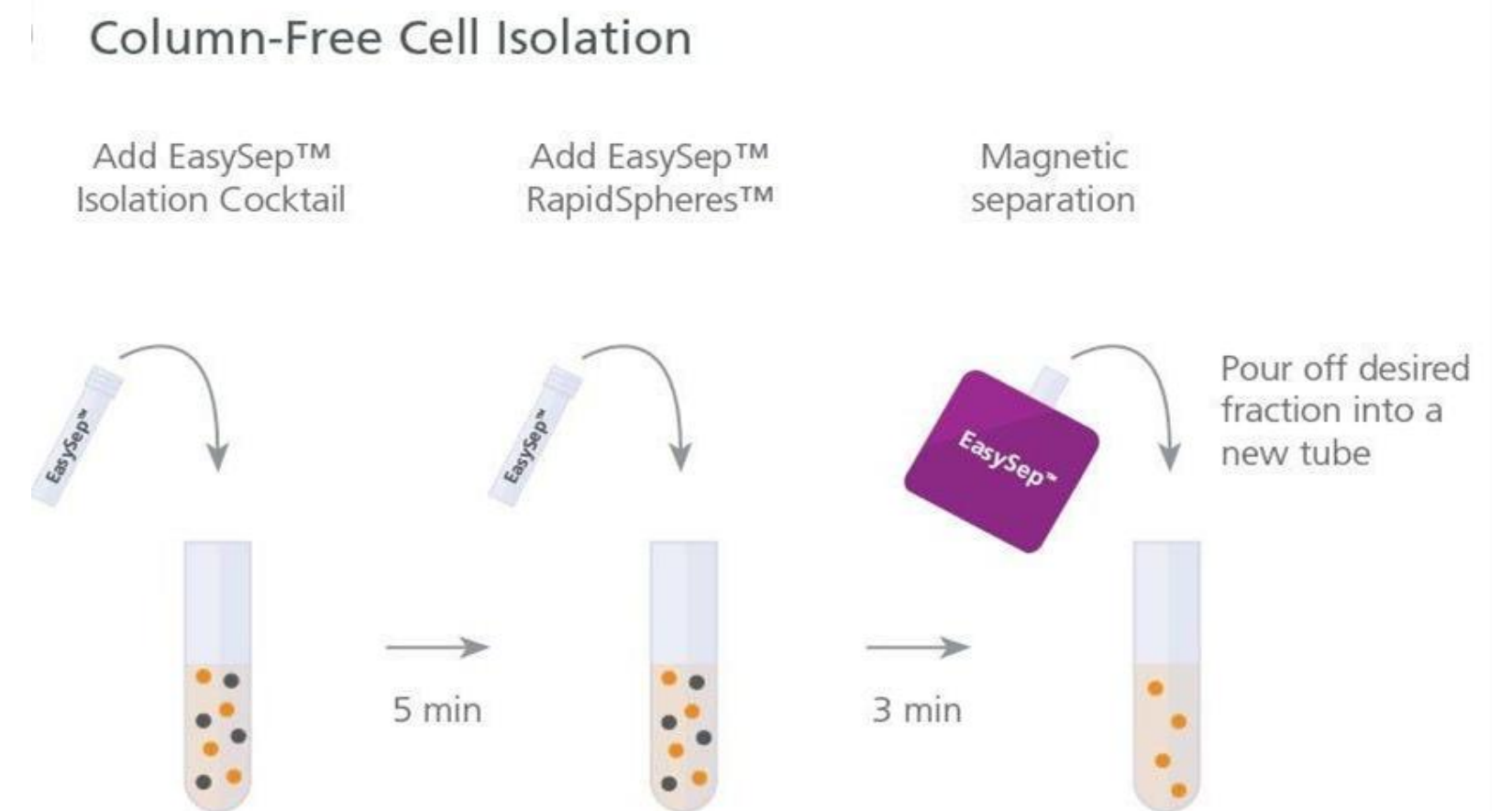

EasySep™ Product Portfolio Overview

Human Cell Separation Kits

- Isolation, Positive & Direct formats
- PBMCs, T/B/NK cells, whole blood
- For immunology & translational research



EasySep™ Human Extracellular Vesicle (CD61) Depletion and Positive Selection Kit

Mouse & Rat Cell Separation Kits

- Isolation & Positive selection formats
- Hematopoietic progenitors, Tregs
- For immunology & stem cell models



EasySep™ Mouse T Cell Isolation Kit

Non-Human Primate Kits

- Isolation & Positive selection
- CD4+ T cells, Pan-T cells
- For preclinical & vaccine research



EasySep™ Non-Human Primate T Cell Isolation Kit

Specialized Formats & Technologies

Direct Kits

- Isolate cells directly from whole blood
- No preprocessing required
- Simplifies workflow



EasySep™ Direct Human B Cell Isolation Kit

Release Kits

- Remove magnetic particles post-isolation
- Obtain truly particle-free cells
- Ideal for functional assays



EasySep™ Release Human APC Positive Selection Kit

Automation-Compatible

- Pre-formatted for RoboSep™ instruments
- High-throughput and reproducible
- For standardized workflows



EasySep™ Human CD4+ T Cell Isolation Kit

EasySep™ in In Vitro Diagnostics



SepMate™-15 (IVD) / SepMate™-50 (IVD) tubes

- Density gradient isolation of PBMCs in IVD workflows
- Manufactured under cGMP



EasySep™ Human Bone Marrow CD138 Positive Selection Kit (IVD)

- Enrichment of CD138+ plasma cells from bone marrow for diagnostic assays
- First hematopoietic cell enrichment kit classed as an IVD by FDA
- Used upstream of assays like FISH to improve sensitivity in multiple myeloma diagnostics.
- Approved / registered in US, EU, UK, Canada, etc.



EasySep™ Buffer (IVD)

- Buffer used in IVD-labeled EasySep™ kits to maintain consistency

EasySep™ Magnets



EasySep™ Magnet

- Standard small-volume isolations (up to 2.5 mL)
- Compact, single-tube format, ideal for most kits



EasyEights™ Magnet

- Process up to 8 samples in parallel
- Increases throughput and consistency



“The Big Easy” EasySep™ Magnet

- Large-volume samples (up to 14 mL)
- Handles bigger sample sizes and no columns needed



Easy 250 EasySep™ Magnet


- Very large-scale processing (up to 250 mL)
- For high-volume or clinical research settings

Key Advantages of EasySep™


- **Speed:** < 20 min isolation time
- **Purity & Viability:** > 94% purity, gentle on cells
- **Scalability:** from small research to clinical-scale kits
- **Ease of Use:** minimal training needed
- Compatible with FACS: Cell Enrichment
- **Automation Ready:** compatible with **RoboSep™** instrument

Downstreaming Processes

- Functional cell assays eg. activation and cytotoxicity assays
- Flowcytometry
- Cell culture and expansion
- Molecular and Omics analysis
- Cell therapy and diagnostics
- Cryopreservation

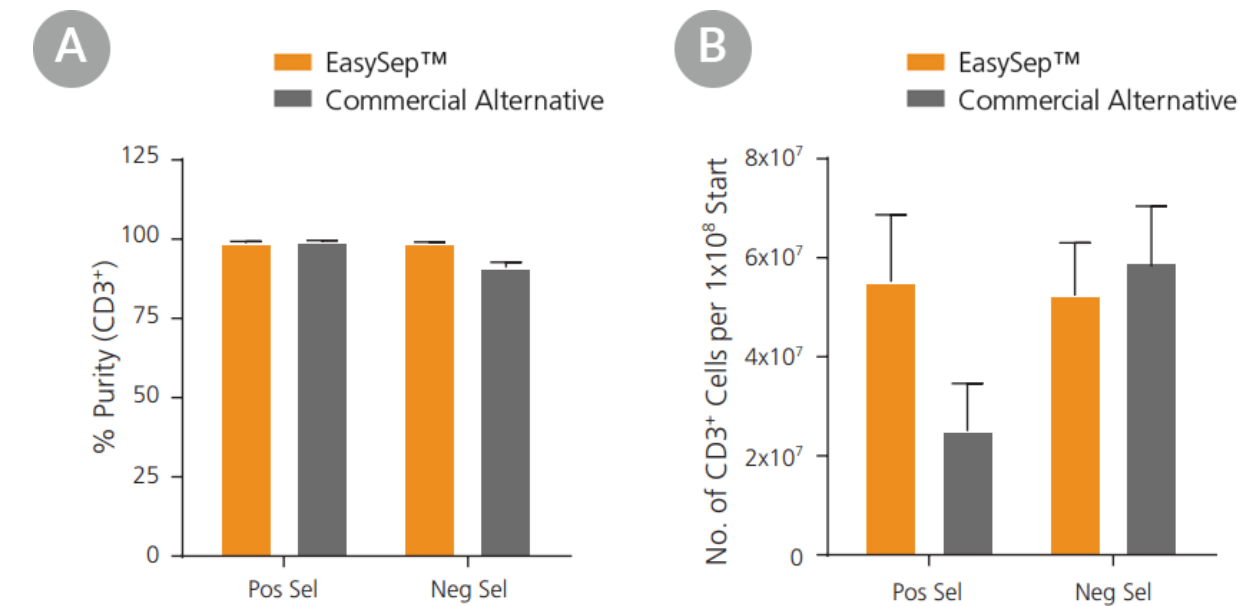

Purity and Recovery of EasySep™ Isolated T Cells

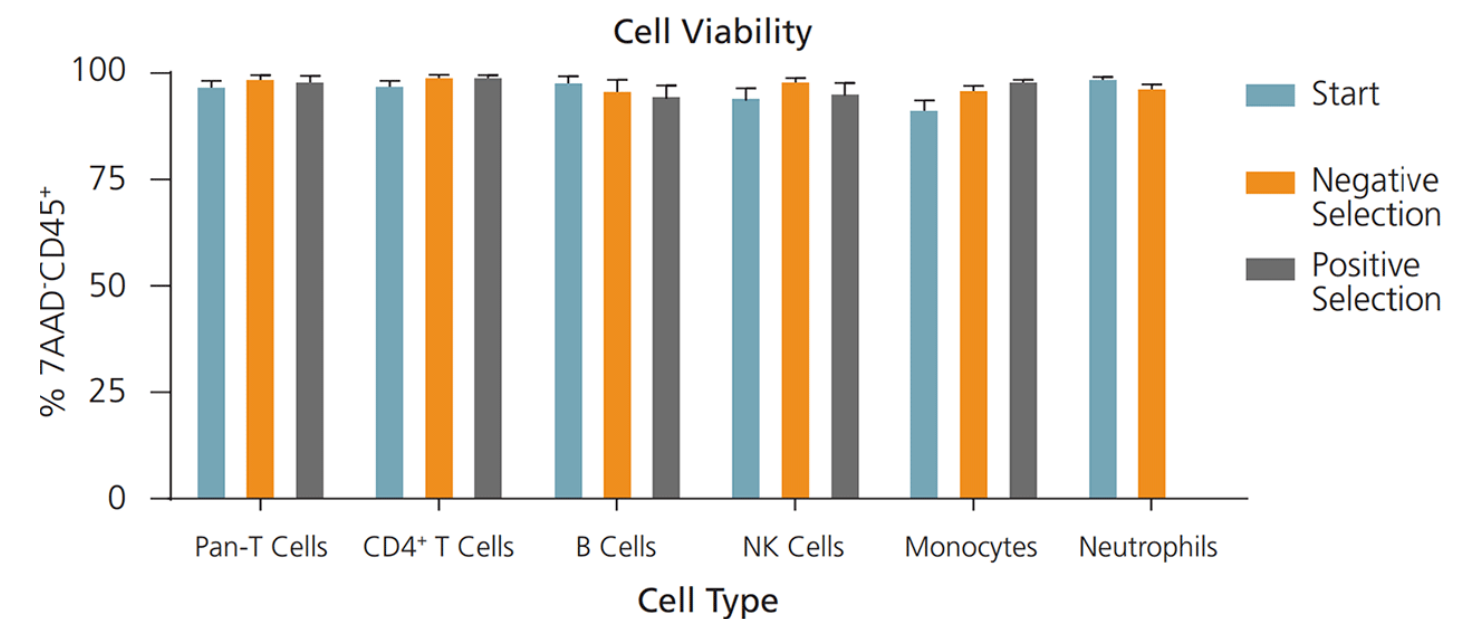


Viability of EasySep™ Isolated Cells

How STEMCELL Supports Customer Needs & Partners in Research



Support

Expert PhD-level technical assistance

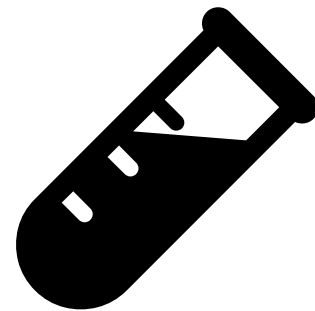


Quality

ISO- and GMP-compliant manufacturing



Workflow

From cell isolation (EasySep™) to culture (MesenCult™, TeSR™) to analysis

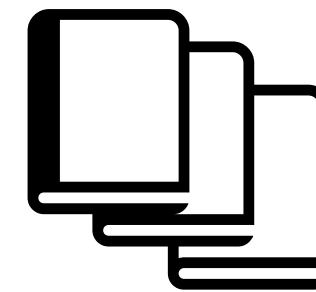

Traning

Webinars, protocols, and workshops




Collaboration

long-term partnerships to accelerate translational research.

How I'd Leverage the Portfolio in Outreach

Personalizing Solutions to Address Client Needs

- Ask prospect: "What cell type + species do you work with?" → tailor kit suggestion
- Use portfolio breadth as differentiation: "We have a kit for that too."
- Cross-sell accessories / magnet upgrades / RoboSep versions
- Educate on use-case comparisons (release vs negative vs direct)

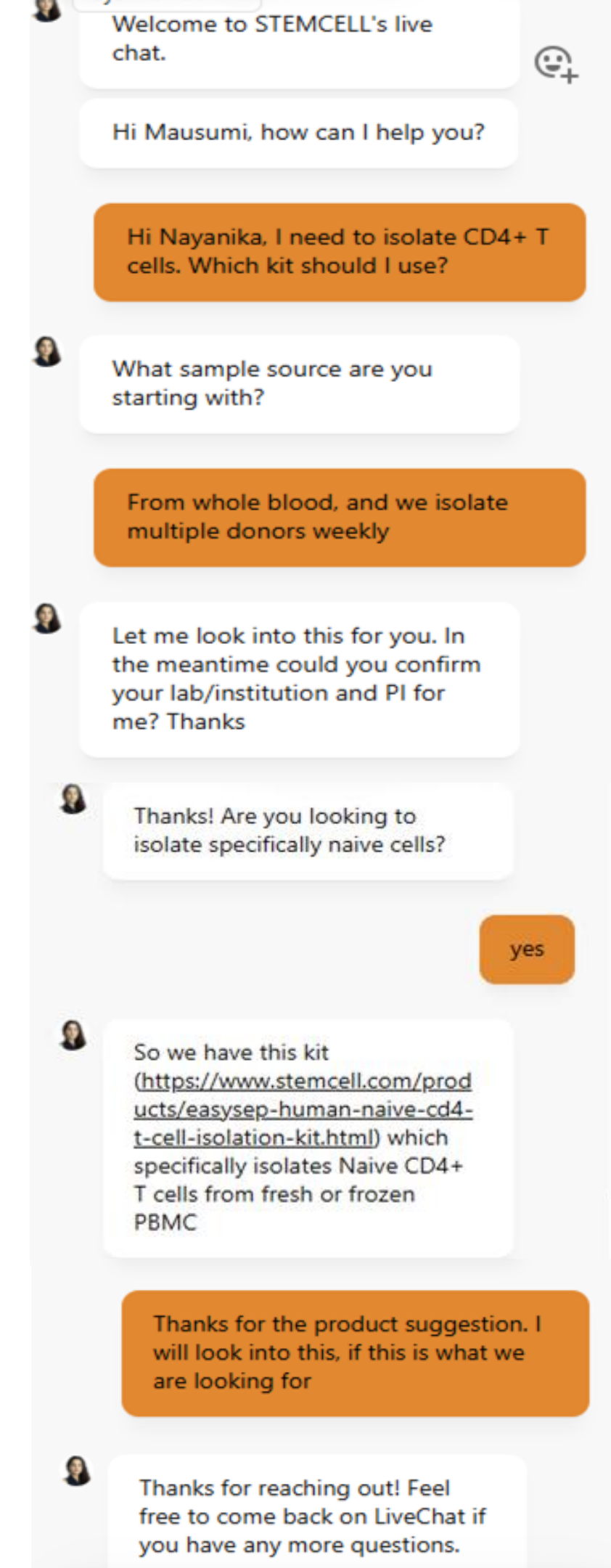

Acknowledgment



All images taken from STEMCELL Technologies website

Thank You for Your Attention!